ADHESIX[®] BIORING[®]

ADJUSTABLE GASTRIC BANDING & ADHESIVE PORT

en	Instructions for use	Page	2
fr	Notice d'instructions	Page	4
de	Gebrauchsanweisung	Seite	6
it	Istruzioni per l'uso	Pagina	8
es	Instrucciones de uso	Pagina	10
ελ	Οδηγίες χρήσης	σελίδα	12
hu	Használati utasítás	Összefoglalás	14

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Version du 27/06/2018

ADHESIX[®] BIORING[®] DEVICE

Adjustable Gastric Banding & adhesive port

Gastric ring kit + adhesive port

DESCRIPTION

The **ADHESIX[®] BIORING[®]** adjustable gastric band is composed of:

- A radiopaque body, comprising an external reinforcement and an inflatable inner part. The extremities of the body form the closure system of the band after its placing around the stomach.
- A radiopaque injection port, which is usually placed in the left subcostal or presternal region of the patient. This site enables the percutaneous adjustment of the diameter of the adjustable gastric band.
- A radiopaque catheter, which links the body of the adjustable gastric band to the injection port. The ensemble is non-absorbable.
- A stiffener, which secures the connection between the injection port and the catheter via a click-shut system.

COMPOSITION

The body is composed of dimethylsiloxane. The external reinforcement is made of polyester which is fully impregnated with dimethylsiloxane, barium sulphate, and a biocompatible dye.

The catheter is made of dimethylsiloxane impregnated with barium sulphate, and printed with biocompatible ink.

The injection port is composed of polyoxymethylene impregnated with barium sulphate. Its septum is composed of dimethylsiloxane. A textile part is fixed to the port by means of a thread in polyethylene terephthalate. This knitting is in polypropylene soaked with an adhesive resorbable hydrogel (Polyvinylpyrrolidone – Polyethylene glycol).

The stiffener is made of dimethylsiloxane and barium sulphate.

No human or animal origin. Semi-resorbable.

INDICATIONS

The **ADHESIX[®] BIORING[®]** adjustable gastric band system is indicated for the reduction of excess body weight in persons affected with morbid obesity, and who are resistant to dieting, i.e. persons whose Body Mass Index (BMI, formula = weight in kilograms / (height in metres)²) is greater than or equal to 40 kg/m², or 35 to 40 kg/m² when combined with at least one comorbidity.

The surgeon is free to prescribe the adjustable gastric band at his own discretion.

CONTRAINDICATIONS

- Inflammatory disease or congenital or acquired anomalies of the gastrointestinal tract.
- Severe organic disease.
- Potential bleeding disorder of the upper gastrointestinal tract or patients who regularly or systematically take aspirin or non-steroidal anti-inflammatory drugs.
- Diagnosis or pre-existing symptoms or family history of connective tissue autoimmune disease.
- Pregnancy.
- Patients with a drug or alcohol addiction.
- Patients who refuse to accept the dietary restrictions required by the procedure.
- Serious psychiatric disorders.
- Do not implant the device at the site of infection or in a septic environment.
- Do not implant the device in the event of allergy to the components.
- Do not implant the device in growing children.

POSSIBLE COMPLICATIONS

Complications that may arise from the use of an adjustable gastric band include risks associated with the anaesthesia, the surgical technique, and the degree of intolerance of the patient to any foreign body implanted in the body.

Perioperative complications (due to the surgical technique):

- Lesions of the spleen
- Perforation of the stomach or oesophagus
- Other bleeding

In the immediate postoperative period:

- | | |
|---|--|
| <ul style="list-style-type: none"> - Nausea, vomiting (especially in the first few days after the intervention, and when the patient eats more than recommended). - Bloating due to gas | <ul style="list-style-type: none"> - Dysphagia, dehydration - Thromboembolic complications - Infectious complications - Postoperative bleeding |
|---|--|

Delayed postoperative complications:

- | | |
|--|---|
| <ul style="list-style-type: none"> - Band dilation - Band erosion - Oesophagal dilation - Pouch dilation - Band alteration - Band intolerance - Pain located at the site of port implantation - Band infection | <ul style="list-style-type: none"> - Slipping of the band, - Gastroesophageal reflux, - Pyrosis, - Gastritis of the upper pouch, - Migration of the band, - Displacement of the injection port, - Infection of the injection port, - Constipation |
|--|---|
- Constrictive obstruction (due for example to oedema, food, an inadequate calibration, slippage of the band, or to torsion or dilation of the upper pouch),
 - Deflation of the band (due to a leak in the band, injection port, or connection tube),
 - Rapid weight loss may lead to symptoms of malnutrition, anaemia, and associated complications

WARNING

The laparoscopic implantation technique is a state-of-the-art procedure. Surgeons who intend to perform this procedure should possess the necessary equipment and have extensive experience in laparoscopy.

PRECAUTIONS

- Store in a dry place, away from light and at room temperature in its original packaging.
- Check the integrity of the packaging (do not use if the labels and / or protective packaging are damaged).
- Do not use if the device is damaged or has past its expiry date.

DO NOT REUSE

DO NOT RE-STERILIZE

SURGICAL TECHNIQUE

Material provided

- One adjustable, radiopaque gastric band, with closure system
 - One syringe
 - One Huber needle
 - One adhesive radiopaque injection port
 - One radiopaque catheter flushing device
 - One handling tag
 - One radiopaque stiffener
- A gastric tube is available separately for the calibration of the small stomach

Preparation of the material: *All of the procedures should be performed in a sterile field.*

- Flush the band with physiological saline using the syringe and catheter rinsing device provided to completely evacuate the air from the device.
- Inspect the inflatable part of the band to detect any irregular inflation or leaks. When handling the band, make sure that it does not come into contact with sharp or pointed instruments.
- Clamp the catheter to prevent air from entering the device again.
- Place the handling tag at the free end of the catheter.
- Flush the implantable injection port with physiological saline, using the syringe and Huber needle.

Furthermore, to get a maximal internal diameter of the band, it is advisable to aspirate the content of the adjustable gastric band completely, using the syringe and catheter-rinsing device, especially during procedures using the "pars flaccida" technique, and in patients with a significant quantity of perigastric adipose tissue.


Surgical technique:

Prophylactic Antibiotics: It is recommended to deliver an antibiotic prophylaxis against skin and gut flora, prior to beginning the procedure.

Various surgical techniques may be used, including: the perigastric technique, the "pars flaccida" technique, or a combination of these two techniques.

Adjustable gastric bands are usually placed laparoscopically, under general anesthesia:

- The anaesthetist introduces the cuffed gastric tube into the stomach, making it possible to calibrate the small stomach.
- Dissect the lesser curvature or the right crus.
- 2nd window next to the His angle.
- Retro-gastric tunnelling.
- Introduction of the adjustable gastric band into the abdomen. When passed through the trocar (12mm minimum), the clip must be positioned on the reinforced part of the buckle.
- Inflate the cuff of the gastric tube with 15 cm³ (ml) of physiological saline.
- Correctly position the band.
- Suture the band, using 3 to 6 sutures, to the gastric serosa, upstream and downstream (anterior gastro gastric suture enclosing the band in a "tunnel").
- Connection of the tube to the injection port.
- General port placement recommendations:
 - Make sure the catheter shows no sharp bending
 - When connecting the catheter to the port make sure you don't create twisting/spring effect on the catheter: check that the port remains totally flat on a surface
 - Make sure the stiffener is properly clipped to the port (click-shut system).
 - In the abdominal wall layer, the path of the catheter must be oblique rather than perpendicular to the port. This prevents lifting effect on the port which with time could lead to port rotation:

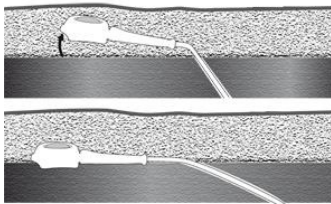Wrong catheter pathway: the catheter may generate a lifting force on the port and lead to port rotation.

Correct catheter pathway: the catheter is in an oblique path through the muscle and generates no lifting force on the port which remains flat on the muscle.

The injection port is self-adhesive. It is not necessary to use sutures; simply position it on the muscle:

- inserted in the aponeurosis of the rectus abdominis muscle under the left costal edge,
- inserted subcutaneously under the left costal edge,
- inserted pre-sternally, on the xiphoid apophysis;
- in any other position according to the surgeon's appreciation.

The bands may also be implanted via laparotomy.

The band can be inflated immediately using physiological saline (NaCl 0.9%) The volume, as well as the interval before filling, is at the discretion of the practitioner/surgeon (maximum volume, see the next table).

Inflating of the band: maximum volume of physiological saline according to sizes of the bands

Band sizes (diameter : mm)	Maximum volume of physiological saline (cm³)
Ø 21/24	10
Ø 27/29	10
Ø 36/38	15

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