BIOMESH® P1

STERILE EXTRAPERITONEAL NON-RESORBABLE PARIETAL REINFORCEMENT IMPLANTS

en	Instructions for use	Page	2
fr	Notices d'instructions	Page	7
de	Gebrauchsanweisung	Seite	13
it	Istruzioni per l'uso	Página	19
es	Instrucciones de uso	Página	25
pt	Instruções de utilisação	Pagina	31
ελ	Οδηγίες χρήσης	σελίδα	37
bg	Инструкции за употреба	Страница	43
da	Brugsvejledning	Side	49
et	Kasutusjuhend	lehel	54
fi	Käyttöohjeet	Sivu	59
hu	Felhasználási Útmutató	Oldal	65
lv	Lietošanas instrukcija	Lappuse	71
nl	Gebruiksinstructies	Pagina	77
ro	Instrucţiuni de Utilizare	Pagină	83
tr	Kullanim Kilavuzu	Sayfa	89
sv	Bruksanvisning	Sida	94
cs	Návod k použití	Strana	99

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CE 1639

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NOT340_260828
Version du 28/08/2026

This release is the last update of the instructions of use and replace the previous edition

Product under regulation 2017/745 (EU).

CE PENDING

Instructions for use *en*

BIOMESH® P1
STERILE EXTRAPERITONEAL NON-RESORBABLE PARIETAL REINFORCEMENT IMPLANTS
SINGLE USE STERILE PRODUCT

MANUFACTURER IDENTIFICATION

COUSIN BIOTECH
ALLEE DES ROSES
59117 WERVICQ SUD
FRANCE

IDENTIFICATION AND DESCRIPTION OF THE DEVICE

BIOMESH® P1 implants are non-absorbable parietal reinforcement implants with:

- Shape memory
- Resistance to suturing
- Flexibility

The blue marking on preshaped (FBIOP13DXX) meshes is the anatomic landmark for the pubic bone and medial side of the patient.

INDICATIONS

Repair and parietal reinforcement of inguinal, femoral, and ventral hernias.

CONTRAINDICATIONS

Do not use in the following cases:

- Allergy to any of the components
- Infected site
- Pregnancy
- Growing children
- Anticoagulant therapy

UNDESIRABLE SIDE EFFECTS

Like any implantable medical device, this implant is susceptible to generate possible undesirable effects:

- Discomfort/Pain	- Adhesion formation	- Hematoma	- Foreign Body Reaction
- Infection	- Fistula formation	- Mesh migration	- Irritation near by organ
- Recurrence	- Erosion	- Mesh deformation	- Ileus / Obstruction
- Inflammation	- Seroma	- Allergic reaction	- Lymphocele
- Urinary retention	- Organ injury	- Chronic Pain	

TARGET POPULATION

- Adult male or female suffering from groin hernias with an indication of surgical repair with extraperitoneal implantation. These include:
 - Inguinal hernias
 - Femoral hernias
- Adult male or female suffering from a symptomatic ventral hernia necessitating a surgical repair with a method that includes an extraperitoneal implantation.

IMPLANTED MATERIALS

Polypropylene
Blue landmark (colorant [phthalocyaninato(2-)] copper) polypropylene marking depending on references (See on the package).

CE PENDING

References	Weight of Polypropylene (**g**)* Theoretical calculation	Weight of blue landmark (colorant [phthalocyaninato(2-)] copper) (**mg**)*
FBIOP10510	0.50	Not applicable for this reference
FBIOP10611	0.66	Not applicable for this reference
FBIOP10715	1.05	Not applicable for this reference
FBIOP11015	1.50	Not applicable for this reference
FBIOP11515	2.25	Not applicable for this reference
FBIOP12536	9.00	Not applicable for this reference
FBIOP11530	4.50	Not applicable for this reference
FBIOP13030	9.00	Not applicable for this reference
FBIOP14545	20.25	Not applicable for this reference
FBIOP16011	0.60	Not applicable for this reference
FBIOP16135	0.74	Not applicable for this reference
FBIOP17515	0.95	Not applicable for this reference
FBIOP1T510	0.51	Not applicable for this reference
XBIOP11515	2.51 per mesh (box of 10 meshes)	Not applicable for this reference
VBIOP13030	9.90 per mesh (box of 5 meshes)	Not applicable for this reference
FBIOP1PR51	0.76	Not applicable for this reference
VCBP11515X	2.51 per mesh (box of 10 meshes)	Not applicable for this reference
VCBP16011X	0.51 per mesh (box of 10 meshes)	Not applicable for this reference
FBIOP13DSL	1.23	≤3.60
FBIOP13DSR	1.23	≤3.60
FBIOP13DML	1.45	≤3.60
FBIOP13DMR	1.45	≤3.60
FBIOP13DLL	1.72	≤3.60
FBIOP13DLR	1.72	≤3.60

* These values are subject to +/- 15% variations

Origin neither human nor animal – Non resorbable.

PERFORMANCES

BIOMESH® P1 allows parietal reinforcement and repair of parietal hernias.

BIOMESH® P1 is particularly adapted to laparoscopic surgery, and to conventional open surgery. This range is biocompatible and designed for extraperitoneal implantation only.

LINK TO SUMMARY OF SAFETY AND PERFORMANCE CHARACTERISTICS:

https://www.cousin-biotech.com/d-6318526d8b749_under-construction.pdf

MRI SAFETY

The implant has not been tested for safety and compatibility in the case of postoperative investigation using Magnetic Resonance Imaging (MRI). However, the implant is mainly made with polypropylene and non-metallic component. Consequently, **BIOMESH® P1** is considered MR safe.

PRECAUTIONS FOR USE

BIOMESH® P1 meshes are delivered sterile (ethylene oxide sterilization).

Before any use, inspect the integrity of packaging and device (of which peelable pouches). Do not use if the device is out of date.

Do not use in the event of deterioration of the device and/or the packaging.

Unpacking and manipulation of the implant must be carried out following asepsis standards in order to guarantee the sterility.

When opening the carton box, the surgeon has to make sure that the lot number, the model and the size of the mesh are equal to the data on the label of the inner pouch. Otherwise, the use of the prosthesis is strictly prohibited.

Prevent the mesh from any contact with objects that may alter their surface.

Surgical use: Temperature from +32°C up to +42°C. Implants should be resistant to biological fluids and to surrounding anatomical tissues. Therefore, the **BIOMESH® P1** has been validated for a long-term implantation.

COUSIN BIOTECH does not offer any guarantee or recommendation as far as the use of a particular type of fixation.

When the surgeon considers it necessary to do so, it is recommended to cut the implant outside the surgical site to avoid any risk of releasing particles into the patient.

This device must be implanted only by a qualified surgeon trained on the use of the product (knowledge of anatomy and visceral surgery).

CE PENDING

IMPORTANT: DO NOT REUSE – DO NOT RESTERILIZE

As specified on the product's labelling, the **BIOMESH® P1** mesh is for single use only. It cannot be re-used and/or resterilized (potential risks would be and are not limited to: loss of the product's sterility, risk of infection, loss of the product's efficiency, recurrence).

INSTRUCTIONS TO BE GIVEN TO THE PATIENTS BY THE SURGEON

The surgeon should inform the patient of the potential physical and psychological restrictions and consequences of implanting the device. The patient must be informed of the surgical risks and possible side effects. The surgeon should invite the patient to return for further consultation if the patient presents symptoms that appear abnormal. An implant card for the patient is provided by the healthcare professional at the clinic. This implant card provides information, for the patient, to identify the device and traceability elements as well as the name, address and website of the manufacturer. The surgeon shall invite the patient to scan the implant card immediately upon receipt to keep track of it in case of loss. The surgeon shall also inform the patient that the summary of the safety and performance characteristics of the device and the instructions for use can be found on the COUSIN BIOTECH website.

STORAGE PRECAUTIONS

Store in a dry place away from sunlight and at room temperature in its original packaging.

SURGICAL TECHNIQUE

Open inguinal hernia repair (Lichtenstein)

➢ The patient is placed in the supine position
➢ Local or general anaesthesia
➢ Iliac incision of approximately 5 cm
➢ Dissection is performed accordingly to the Lichtenstein technique
➢ The hernia sac is identified and reduced
➢ The transversalis fascia can be tightened with a non-absorbable suture. With direct hernias it helps to avoid any immediate post-operative recurrence
➢ The mesh can be cut to fit the patient's anatomy
➢ The prosthesis is progressively put in place
o begin by placing the round tip of the mesh towards the pubic bone
o the upper part of the prosthesis is placed on the conjoint tendon
o carry placing the prosthesis towards the deep inguinal ring so that the first arm is placed behind
o then the second arm is placed so that it encircles the spermatic cord
o Fix the mesh according to surgeon's preference
▪ apex is fixed to the pubic tubercle
▪ the lower border of the mesh is fixed to the free edge of the inguinal ligament
▪ continuous suture extends up just medial to the anterior superior iliac spine
▪ fix the two tails of the mesh together around the spermatic cord
▪ the infero-medial corner of the mesh is attached overlapping the pubic tubercle.
▪ the mesh is fixed to the conjoined tendon
o the cord then rests on the prosthesis
o suction drain can be placed beneath the external oblique aponeurosis, especially in large inguinal hernias, where an extensive dissection was performed
o the aponeurosis of external oblique is then closed with absorbable sutures
o before closure of the surgical incision, its edges are infiltrated with a long-acting local anaesthetic, such as Naropein®
o During peri-operative care of the patient, prophylactic antibiosis is usually given for 48 – 72 hours postoperatively

Laparoscopic inguinal repair (TAPP or TEP)

1/ Example of a TAPP repair
➢ General anaesthesia
➢ A pneumoperitoneum is created in the usual fashion
➢ intra-abdominal findings are reported [intra-abdominal pathology and inguinal hernia defects and sacs].
➢ The two additional trocars are inserted under direct vision.
➢ Dissection is initiated and carried out following classical TAPP technique
➢ The indirect inguinal hernia sac should be dissected carefully from the Spermatic Cord. Direct hernia sacs are easily dissected.

PENDING

➢ The mesh is inserted through the Trocar into the intra-abdominal cavity and deployed over the inguinal region. The blue landmark on the preshaped meshes is the anatomic landmark for the pubic bone and Medial side of the patient. The mesh is attached to Cooper's Ligament, around and lateral to the Inferior Epigastric Vessels using fixations according to the Surgeon's usual technique. The operator should be meticulous to avoid the triangle of doom, triangle of pain, iliac vessels and to place potential fixations lateral to the inguinal ring. Mesh fixation should be at least 1cm from the edge of the mesh with 1 cm spacing between fixation points.
➢ The operator should check the mesh is well anchored to the surrounding structures.
➢ The peritoneum is closed meticulously and no defect between the peritoneum and the abdominal wall should be left open. In addition, it should cover the entire Mesh.
➢ The trocars are removed under direct vision. The fascia of the sub-umbilical trocar site is closed as needed.

2/ Example of a TEP repair

➢ General anaesthesia
➢ Create the pro-peritoneal space.
➢ The 10 mm trocar is inserted.
➢ Create the pneumo-pro-peritoneum.
➢ Two 5 mm ports are inserted under direct vision.
➢ If a peritoneal defect is created at this time it should be immediately closed.
➢ The anatomical landmarks are identical to when performing a TAPP repair.
➢ The indirect inguinal hernia sac should be dissected carefully from the Spermatic Cord. Direct hernia sacs are easily dissected. This is done by gentle traction with atraumatic graspers.
➢ The mesh is inserted into the pre-peritoneal space and deployed over the inguinal region. The blue landmark on the preshaped meshes is the anatomic landmark for the pubic bone and Medial side of the patient.
➢ The mesh is first secured to Cooper's Ligament, and the superior aspect of the pubic ramus. The mesh is then fixed on the posterior aspect of the Linea Alba. The anchoring is continued around and lateral to the Inferior Epigastric Vessels. A few fixations are used to fix the mesh lateral to the internal ring. Mesh fixation points should be at least 1cm from the edge of the mesh with 1 cm spacing between fixation points.
➢ The mesh should gently and generously cover the Iliac Vessels without major gaps.
➢ The operator should check the Mesh is well anchored to the surrounding structures. The Mesh should not migrate and remain in place.
➢ The pre-peritoneal space will be checked for any peritoneal defects. If any, they should be closed. The 5 mm trocars are removed under direct vision to ensure there is no bleeding. The mesh must lay flat
➢ If there are any questions about a missed peritoneal defect, a completion laparoscopy should be performed.
➢ Trocars are removed and the fascial defect closed with the appropriate suture. The skin edges are approximated in the usual manner.

Extra-Peritoneal Ventral Hernia repair

Example: Retrorectus mesh repair of a ventral hernia (Rives-Stoppa technique)

This technique utilizes the hernia sac to separate the mesh from the intra-abdominal contents.
- Superior to the umbilicus, dissection is performed above the posterior rectus fascia and under the rectus muscle.
- Below the umbilicus, dissection occurs in the preperitoneal space due to the lack of a posterior rectus sheath.
- The mesh is placed in the newly formed pre-peritoneal space.
- Fixation of the mesh if needed, to the muscle layer above
- The peritoneum must be meticulously checked and no defect in the peritoneum should be left open to prevent any contact between the mesh and intra-abdominal contents.
- Redon's drains under aspiration are usually placed. The anterior sheath is closed. If tension is present, relaxing incisions may be used
- Skin closure.
- Using the Redon's drains, fill the pre-peritoneal space with an antiseptic solution and clamp the drains for 1 hour.
- A belt can be used to wrap the patient's abdomen for a few days to reinforce the abdominal wall.

EXPLANTATION AND DISPOSAL OF DEVICES

This device is intended for permanent parietal reinforcement and to remain permanently implanted. Routine removal or replacement is not required. In the event of certain complications, a secondary surgical procedure for partial or complete mesh explantation may be considered, if clinically indicated. Complete removal of all implant material may

CE PENDING

not always be technically feasible. The decision to explant the device lies with the responsible surgeon and shall be based on the individual clinical situation and current medical practice.
Mesh explantation is a technically demanding reoperation and should only be performed by a surgeon who is appropriately qualified, trained, and experienced in mesh removal techniques. As with any reoperation, mesh explantation carries significant risks, including but not limited to bleeding, infection, injury to adjacent anatomical structures (e.g. nerves, vessels, bowel, bladder, or other organs), anaesthesia related complications, recurrence or persistence of the original condition, worsening of symptoms, chronic pain, and the need for additional surgical interventions. Patients should be informed that mesh explantation may resolve, improve, partially improve, fail to resolve or worsen the reported symptoms or adverse effects.
If a hernia mesh implant is explanted and retrieval data are to be collected, the device should be retrieved and handled in accordance with ISO 12891 1:2015 “Implants for surgery – Retrieval and analysis of surgical implants – Part 1: Retrieval and handling”. Explanted devices that are not returned for analysis must be disposed of in accordance with applicable local regulations for infectious and hazardous clinical waste.

IMPORTANT

For more information about the use of this product, please contact your COUSIN BIOTECH representative or distributor.

INFORMATION REQUEST AND COMPLAINTS

In accordance with its quality policy, COUSIN BIOTECH is committed to produce and supply a high-quality medical device. However, if health professionals (client, user, prescriber, etc.) has a complaint or cause for dissatisfaction with product in terms of quality, safety or performance, they must inform COUSIN BIOTECH as soon as possible. In the event of implant failure or if the implant is implicated in serious adverse effects experienced by the patient, the health centre must follow the legal procedures in the country in question, and inform COUSIN BIOTECH without delay. For any correspondence, please specifiy the reference, batch number, reference contact details and a comprehensive description of the incident or complaint. If further information is needed or required, please contact your COUSIN BIOTECH representative or distributor.

Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.


IMPLANT CARD INFORMATION LEAFLET

How to complete the implant card ?

The implant card, for the patient, should be completed as follows:

Front side of the implant card

The user has to complete the implant card before transmitting the implant card to the patient. Please complete it as follows:

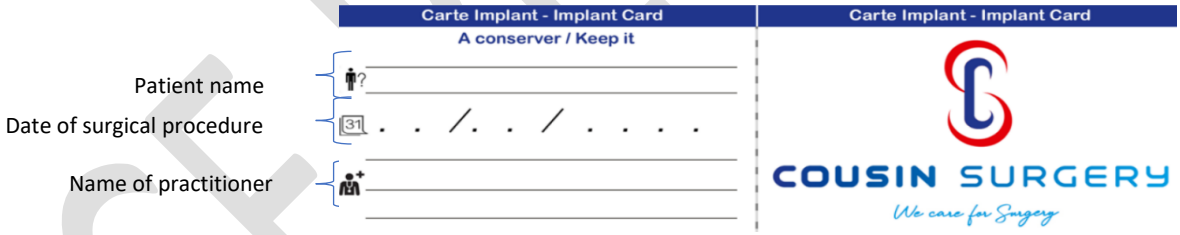Back side of the implant card

No fields need to be completed by the user. The following information on the device is printed directly by the manufacturer on the back of the implant card. A QR code allowing the patient to consult the instruction for use, as well as the website address, are available on the implant card.

Information to identify the device	Traceability information
• Indication that the product is a medical device • Description of the device family • Trade name of the device • Reference of the device • Device is MR safe • Denomination name of the device in sales countries	• UDI number • Batch number • Name and address of the legal manufacturer

Explanation of the used symbols in the Implant Card

All symbols are commonly used and explained at the end of the instruction for use.

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